




ARMOURFIT™

FIT INSTRUCTIONS

ARMOURFIT™ STRAPPED
ARMOURFIT™ STRAPLESS
UA HOOPS
UA SOCCER

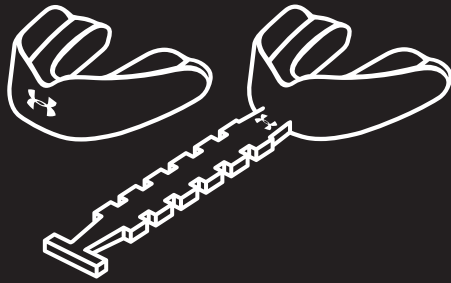

CARE INSTRUCTIONS

- UA Mouthwear Freshener spray is recommended after each use to keep your mouthguard fresh and clean.
- Protect your mouthguard by storing it in the UA Mouthwear case.



Your strapped mouthguard can be easily converted to a Trim-Tab.

1. Fit your mouthguard according to the fit instructions.
2. Securely hold the mouthguard, and using scissors carefully cut off the strap just above the UA logo.

NOTE: Scissors are sharp. Use caution when handling.

STOVE-TOP

READ ALL INSTRUCTIONS PRIOR TO FITTING YOUR MOUTHGUARD.

What you need for a stove-top fitting:

- Pot of boiling water 4" deep
- Slotted spoon
- 9" x 9" pan of cold water
- Mirror
- Timer



1. Bring pot with 4 inches water to a boil and then remove from heat source. **



2. Submerge mouthguard in hot water for 30 seconds. **Strapless** - carefully remove mouthguard from hot water with slotted spoon. **Strapped** - do not submerge the strap.



3. Immediately submerge mouthguard in pan of cold water for no more than 2 seconds. Do not allow sides to fold over as the material will "glue together" making it difficult to get an optimal fit.



4. Gently glide your thumbs inside the molar channel to ensure there is ample space for the mouthguard to fit around your teeth. Using a mirror, carefully place the mouthguard in your mouth, lining up the "notch" in the center of the mouthguard with the center of your front teeth.



5. To ensure proper fit, simultaneously perform below steps for 2 minutes:

- a. Bite down gently and naturally on mouthguard.
- b. Press tongue firmly against the roof of mouth and with light sucking action, create a vacuum chamber by drawing out all water and air.
- c. Press fingers along outside of mouthguard to form to teeth.



6. After fitting process is complete, remove mouthguard from mouth and lay completely flat in 9" x 9" pan of water for 5 minutes to allow the mouthguard to cure/harden.

RE-FITTING

If re-fitting is desired, repeat steps 1-6. If you would like a tighter fit, add 10 seconds to stove-top time in Step #2.

ArmourFit™ is designed to be worn on the upper teeth and must be fitted as instructed to function properly. A damaged, altered or improperly fitted mouthguard will not provide the same protection.

**** CAUTION: BOILING WATER AND ITEMS THAT HAVE BEEN MICROWAVED WILL BE HOT AND CAN CAUSE INJURY.**

Bite Tech • 20 Glover Avenue • Norwalk, CT 06850 • UA approved licensee • Visit www.bitetech.com

MICROWAVE

What you need for a microwave fitting:

- Microwave oven (1100 watts +)
- Wet paper towel
- Microwave-safe dish
- Mirror
- 9" x 9" pan of cold water
- Timer



1. Wrap mouthguard in wet paper towel and place on microwave safe dish.

Strapped - do not wrap the strap.



2. Microwave mouthguard for 45 seconds on high then remove from microwave. **



3. Immediately unwrap mouthguard and submerge in pan of cold water for no more than 2 seconds. Do not allow sides to fold over as the material will "glue together" making it difficult to get an optimal fit.



4. Gently glide your thumbs inside the molar channel to ensure there is ample space for the mouthguard to fit around your teeth. Using a mirror, carefully place the mouthguard in your mouth, lining up the "notch" in the center of the mouthguard with the center of your front teeth.



5. To ensure proper fit, simultaneously perform below steps for 2 minutes:

- a. Bite down gently and naturally on mouthguard.
- b. Press tongue firmly against the roof of mouth and with light sucking action, create a vacuum chamber by drawing out all water and air.
- c. Press fingers along outside of mouthguard to form to teeth.




6. After fitting process is complete, remove mouthguard from mouth and lay completely flat in 9" x 9" pan of water for 5 minutes to allow the mouthguard to cure/harden.

RE-FITTING

If re-fitting is desired, repeat steps 1-6. If you would like a tighter fit, add 15 seconds to microwave time in Step #2.

NEVER MICROWAVE FOR MORE THAN 70 SECONDS

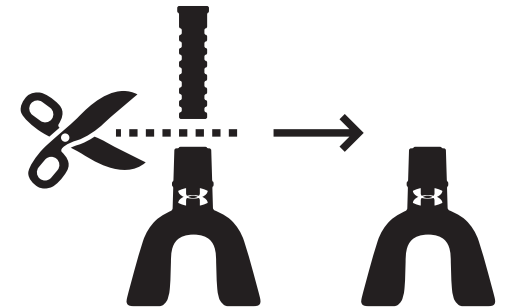

DENTAL LIMITED WARRANTY

CERTIFICATE OF LIMITED WARRANTY AND DISCLAIMER

WARRANTY: Bite Tech, Inc. (Bite Tech) warrants this mouthguard to be free from defects in material and workmanship which could cause damage or injury to the sound, natural teeth of the purchaser herein when properly fitted and used for a period of one (1) year from the date of purchase.

LIMITATION OF LIABILITY: This limited warranty applies only to injuries to the purchaser's teeth that are free from any pre-existing

natural or unnatural defect, including prior dental work. This limited warranty only applies if the purchaser suffered a dental injury in a supervised athletic contest or coach-supervised training session not otherwise excluded. No warranty is provided for any damage or injury incurred during or in connection with any martial arts, mixed martial arts, boxing or similar fighting or combat sport or activity. Bite Tech's obligation under this warranty shall be limited to the actual costs incurred by the user for necessary dental treatment by an appropriately licensed dental professional up to a maximum liability of \$500 per injured or replaced tooth with a maximum liability per incident of \$32,000. Bite Tech's obligation shall be reduced by any insurance coverage available to the user and applicable to the treatment of the injury claim. This limited warranty does not provide payment or indemnity for orthodontic treatments or preparation of a tooth to be used as an anchor for a bridge. The remedies set forth herein are exclusive and the total liability of Bite Tech with respect to use of this mouthguard, whether based on contract, warranty, negligence, indemnification, strict liability or otherwise. In no event shall Bite Tech be liable for consequential, incidental or special damages.

CONDITIONS: This limited warranty applies only to injuries to the purchaser's teeth that are free from any preexisting natural or unnatural defect, including prior dental work. This limited warranty only applies if the purchaser suffered a dental injury in a supervised athletic contest or coach-supervised training session. Bite Tech shall have no responsibility or liability hereunder, unless:

1. the mouthguard was fitted in accordance with the printed instructions that accompany the mouthguard; and
2. the mouthguard was not damaged or altered in any way, whether by chewing on the mouthguard or by other action to alter its configuration or chemistry (other than by the prescribed fitting), except as occurs by natural use; and
3. a complete, written notice of the injury and damage, together with a copy of the receipt of purchase, is received by Bite Tech within thirty (30) days after the occurrence of the injury; and
4. the user provides a written statement from an appropriately licensed dentist providing necessary treatment within thirty (30) days after the date of injury that substantiates the injury, damage and treatment provided, together with identification of insurance coverage and insurance payments; and
5. the user provides a statement by a supervising coach or participating game official that the mouthguard was in proper use at the time of injury; and
6. the applicable mouthguard is returned to Bite Tech with the written notice of injury. **The mouthguard shall become the property of Bite Tech upon receipt of return.**

NO MOUTHGUARD CAN PREVENT ALL INJURIES OR ALL DENTAL DAMAGE AND NO SUCH WARRANTY IS BEING MADE. THIS WARRANTY IS MADE IN LIEU OF ALL OTHER WARRANTIES EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTIES FOR A PARTICULAR PURPOSE AND ANY IMPLIED WARRANTY ARISING OUT OF A COURSE OF DEALING, CUSTOM, USAGE OR TRADE. EXCEPT AS IS OTHERWISE EXPRESSLY PROVIDED HEREIN, BITE TECH DISCLAIMS ALL LIABILITY FOR LOSS, INJURY OR DAMAGE BASED UPON NEGLIGENT DESIGN OR DEFECTS IN MATERIALS AND WORKMANSHIP FOR THE WITHIN MOUTHGUARD. Keep this dental warranty, disclaimer and store receipt in your records. In case of injury, mail mouthguard, store receipt, the notice of injury and accompanying statements to Bite Tech, Inc. This product is protected by one or more U.S. and international patents and pending patent applications. Please visit www.ArmourBite.com for more details.

MAN-01026-001